




replacement parts industries, inc.

CALL 800-221-9723 or 818-882-8611 • FAX 818-882-7028
E-MAIL: techsupport@rpiparts.com • WEBSITE: www.rpiparts.com

RPP819 MOUNTING PAD INSTALLATION INSTRUCTIONS

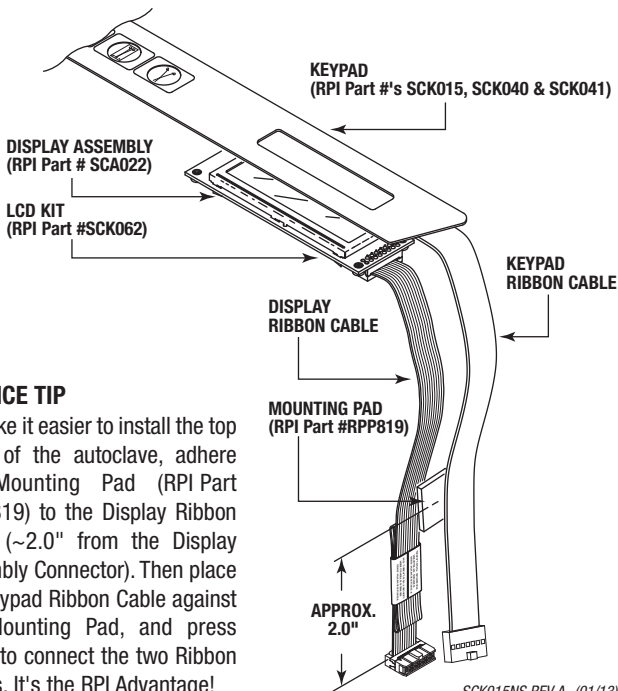

SERVICE TIP

To make it easier to install the top cover of the autoclave, adhere the Mounting Pad (RPI Part #RPP819) to the Display Ribbon Cable (~2.0" from the Display Assembly Connector). Then place the Keypad Ribbon Cable against the Mounting Pad, and press firmly to connect the two Ribbon Cables. It's the RPI Advantage!

SCK015NS REV A- (01/13)